
Quick Start Guide

METINF006 - Non Contact Infrared Thermometer

For the full instructions, visit <u>www.tiptemp.com/METINF006</u>. Click the item, Scroll to 'Support'.

Measure Forehead Temperature

- 1. Install provided (2) AAA batteries
- 2. Make sure probe lens is clean
- 3. Power On: Press 'ON/MEM'
- Measure body temperature on the Forehead: While in Forehead mode (default), hold thermometer within 1.5 inches (4 cm) from central forehead. Press `START' with the `Probe Lens' pointed at the surface. Measurement may take up to 1 second to register.
- 5. Power Off: Manually power off the device by pressing the 'ON/MEM' button. The device will automatically power off if left idle for more than 1 minute.

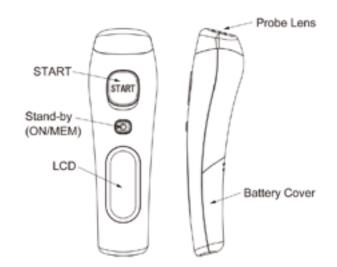

erature Products

NOTE: The thermometer is not meant to replace a visit to the doctor. Please also remember to compare the measurement result to your regular body temperature. Please consult with your doctor if you have health concerns.

Functions		
Function	Details	LCD Icon
Forehead Temperature	*Please see the "Measure Forehead Temperature" section to learn how to measure the body temperature.	*F#
Surface Temperature	Press the 'START' button once to switch to general surface temperature mode. The icon on your display will change.	18.4 °Fm)
Temperature Indicator	If a temperature of 99.5°F (37.5°C) is detected in Forehead Mode, three short beeps followed by one long beep will sound, warning the user of a high temperature.	(00.95)
Memory	There are up to 25 stored temperature measurements. When on, press the 'ON/MEM' button to see the recorded temperature values.	98.25
°F / °C Switch	In 'Power Off' mode, press and hold the 'START' button. Then, press the 'ON/MEM' button for 3 seconds. The '°F' icon will switch to '°C'. NOTE: This will also clear the memory.	2 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0
Mute Mode	When on, press and hold the 'ON/MEM' button for 3 seconds. The icon will flash on the LCD. Release 'ON/MEM' button. NOTE: Holding 'ON/MEM' button for 5 seconds will power off the thermometer without muting.	(of sol

www.tiptemp.com

340 W Broad Street, Burlington, NJ 08016 - 1(800)TIP-TEMP - EveryDegreeMatters@tiptemp.com