


Household Guardian


Infrared Thermometer User Guide



sku SQUEALM003

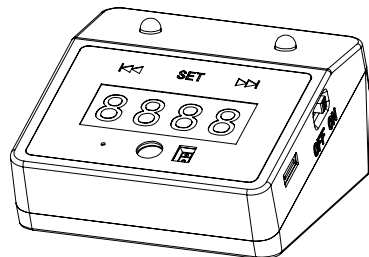

1. Introduction

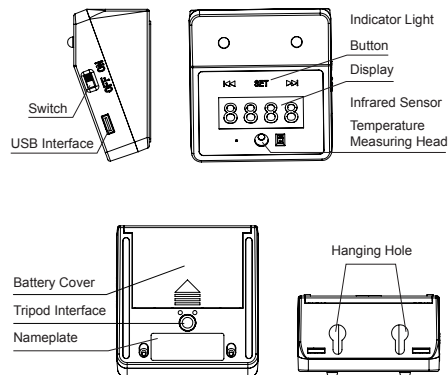

2. Specifications

- Accuracy: $\pm 0.2^{\circ}\text{C}$ ($34\sim 45^{\circ}\text{C}$, place it in the operating environment for 30 minutes before use)
- Abnormal automatic alarm: flashing + "Di Di" sound
- Measuring distance: 3cm~8cm automatic measurement
- Screen: digital display
- Transmission interface: Type-C
- Adaptable power supply: 3 AAA batteries (external power supply 5V1A)
- Install method: Placed on the desktop, nail hook, bracket fixing
- Environment temperature: $10^{\circ}\text{C}\sim 40^{\circ}\text{C}$ (recommended $15^{\circ}\text{C}\sim 35^{\circ}\text{C}$)
- Measuring range: $0\sim 50^{\circ}\text{C}$
- Response time: 0.5s
- Standby: 5S automatic standby during non-working hours
- Input: DC 5V 500mA
- Weight: 73g
- Dimensions: 68x73x36.2mm
- Counting capacity: 9999 times
- Temperature record: 30 sets

3. State Description

- Standby: The display shows the current time.(Fig. 1)
- Insufficient power: The horizontal bar in the display turns on in turn. (Fig. 2)
- Normal temperature: the display shows the temperature and the green light flashes at the top. (Fig. 3)



(Fig. 1)



(Fig. 2)



(Fig. 3)

- Abnormal temperature: the display shows the temperature and the red light flashes at the top.(Fig. 4)



(Fig. 4)

- Out of range:

- 'Lo': Ultra-low temperature than 35 degrees Celsius.(Fig. 5)
- 'Hi': Ultra-high temperature than 42 degrees Celsius.(Fig. 6)



(Fig. 5)




(Fig. 6)

4. Operation Instruction

- Temperature measurement:


Front of the thermometer and approach it with 3-8cm.(Fig. 7)

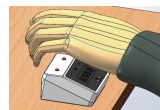


(Fig. 7)

- Temperature unit setting:

Press the SET button to set temperature unit, press the left and right buttons to switch between Celsius and Fahrenheit (Fig. 8, Fig. 9), C: Celsius, F: Fahrenheit.

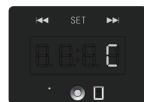

(Fig. 8)

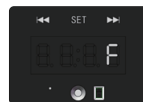


(Fig. 9)

- Mode setting:

Press the SET button twice to set the working mode, press the left and right buttons to switch the mode. Surf :surface mode, body: body mode, coun: counting mode. (Fig. 10, Fig. 11, Fig. 12)



(Fig. 10)

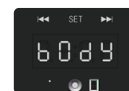


(Fig. 11)



(Fig. 12)


- Alarm temperature setting: Press the SET button three times to enter the alarm temperature setting mode, press the left and right buttons to set it. (Fig. 13)

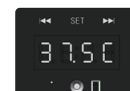

(Fig. 13)

- Time setting:

Press the SET button four times to enter the time setting mode. (Fig. 14)

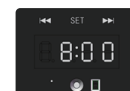

(Fig. 14)

- Battery installation:

- (1) Open the battery cover.
- (2) Install three AAA batteries correctly.



- Install method: Placed on the desktop, nail hook, bracket fixing.



Warnings

- It is the user's responsibility to ensure the electromagnetic compatibility environment of the device so that the device can work normally.
- It is recommended to evaluate the electromagnetic environment before using the device.
- When changing the operating environment, the device must be left for more than 30 minutes.
- Keep away from air conditioners, fans, etc.
- Please use qualified and safety-certified batteries

LimitedWarranty:

1.Warrantor: Dealer, Distributor, Retailer, and Manufacturer

2.Warranties and Remedy

We believe that this is a high quality product. Although we test all products for proper functionality, we cannot guaranty that there will never be a defective unit. For this reason, it must be clear that the Warrantors are not insuring your product/premises or guaranteeing that there will not be damage to your person or property if you use this Product. If this warranty is unacceptable, please return the unused product for a full refund.

One Year Limited Warranty-TIPTEMP warrants this product to be free from defects in material and workmanship under normal use for one year, and is not responsible for consequential damage or installation costs of any nature. In event that the Product does not conform to this Warranty at any time during the period of one year from original purchase date, Warrantor will repair the defect and return it to you at no charge. Important: The Warranty is limited to replacement of the Product ONLY.

This warranty shall terminate and be of no further effect at the time the Product is 1) damaged by extraneous causes such as fire, water, etc. or not maintained as reasonable and necessary: 2) modified: 3)improperly installed: 4) repaired by someone other than the Warrantor: 5) used in a manner or purpose for which the Product was not intended.

WARRANTORS' OBLIGATION UNDER THIS WARRANTY IS LIMITED TO REPAIR OR REPLACEMENT OF THE PRODUCT ONLY. THIS WARRANTY DOES NOT COVER PAYMENT OR PROVIDE FOR THE REIMBURSEMENT OF PAYMENT FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES.

It must be clear that the Warrantors are not insuring your products/premises or guaranteeing that there will not be damage to your person or property if you use this Product. The Warrantors shall not be liable under any circumstances for damage to your person or property or some other person or that person's property by reason of the sale or use of this Product, or its failure to operate in the manner in which it is designed. The Warrantor's liability, if any, shall be limited to the original cost of the Product only. Use of this Product is at your own risk.

3. Procedures for obtaining performance for Warranty:

In the event that the Product does not conform to this Warranty, the Product should be shipped or delivered freight prepaid to a Warrantor with evidence of original purchase. If in any way you are not comfortable with the product or its Limited Warranty, we encourage you to return it unused for a full refund.

DO NOT DISASSEMBLE THIS PRODUCT. It does not contain any user serviceable components.

Attn.: CUSTOMER SERVICE DEPT.

TIP TEMPerature Products

340 W Broad Street

Burlington,NJ 08016

609-239-1900

www.tiptemp.com